



## A.M.I.® RESOFIX™

### Anchoring System for reconstruction of the Anterior Cruciate Ligament!

- Provides highest pull-out forces of Patellar- or Harmstring Ligament Transplants.
- Avoids turning, dislodging or damaging of transplant
- Absorbable
- Put in place fast, easy and safe!

# A.M.I.® RESOFIX™ Anchor

## A.M.I.® RESOFIX™ Anchor

A.M.I.® RESOFIX™ Anchor is a 2-shell expansion bolt used for the fixation of Patellar Ligament- or Harmstring Transplants in ACL-plasties.

The system is fully absorbable and is made of an amorphous Poly DL-lactide resorbable Polyester (R 207 of Boehringer Ingelheim, Germany).

### No damage, no turning, no dislodging of transplant during its fixation!

- High pull out forces = Fast recovery and mobilisation of patient!
- One system for Patellar Ligament- or Harmstring Transplants = Versatile, easy and safe!
- Size of bone-donation site for harvesting of the Patellar Ligament Transplant is reduced by 60 % = Less trauma for the Patella at the donation site!
- Full control of proper placement of transplant before its fixation. = Safe for patient and surgeon.
- Amorphous resorbable material. = No artefacts in post-OP imaging, no removal of metallic-screws or -parts.

#### System / Order Code:

##### Anchors:

RESOFIX™ -S Anchor 10 mm	ARE 201
RESOFIX™ -P Anchor 10 mm	ARE 202

##### Complete Instrument Set


consists of:

- RESOFIX™ Applicator	ARE 200
- RESOFIX™ Dilatation Pin small	ADP 203
- RESOFIX™ Dilatation Pin medium	ADP 209
- RESOFIX™ Dilatation Pin large	ADP 208
- RESOFIX™ Calibration Model	ACM 215
- RESOFIX™ SS-Instrument Tray	AIT 206

#### A.M.I.® RESOFIX™ Anchor

is covered by European and worldwide patents:


EP 0 794 736 in Europe  
US 6.056.750 in the USA

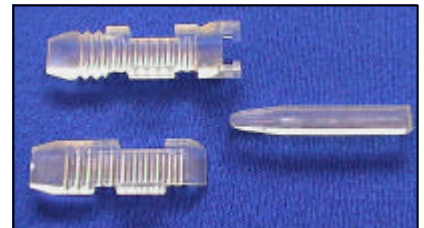


A.M.I.® RESOFIX™ Anchor is a 2-shell expansion bolt.



The fully absorbable A.M.I.® RESOFIX™ Anchor attached to Applicator.

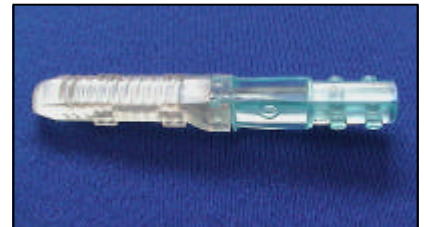

The RESOFIX™ Applicator ARE 200

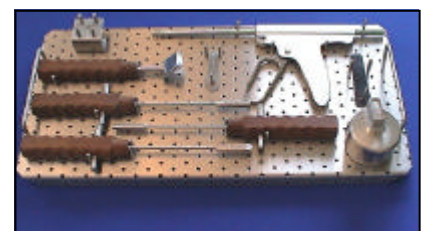

The complete Set of Instruments on the SS-Instrument Tray.

**CJ MEDICAL**

**CJ Medical**

12 Long Wall, Haddenham  
Buckinghamshire HP17 8DL, UK  
Phone and Fax: +44 (1844) 299501  
E-Mail: info@cjmedical.com  
www.cjmedical.com